


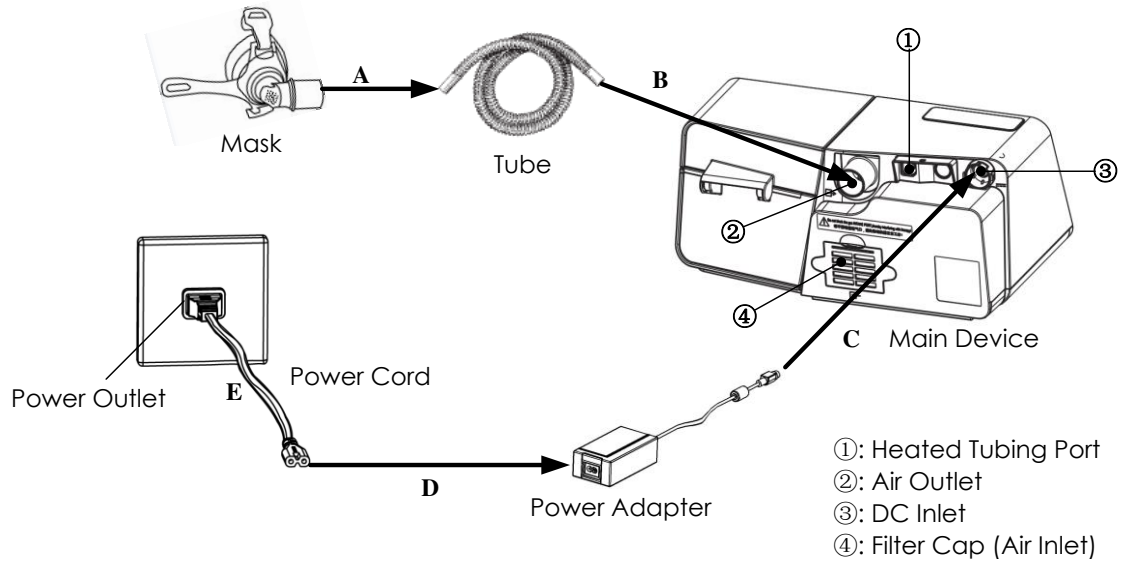


# Luna® G3 BPAP 25A (LG3700) Quick Operation Manual

## 1. Button Functions

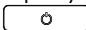
-  Home Button: Return to the previous menu or main interface
-  Knob: Adjust device settings
-  Start / Stop Button: Start / Stop delivering air

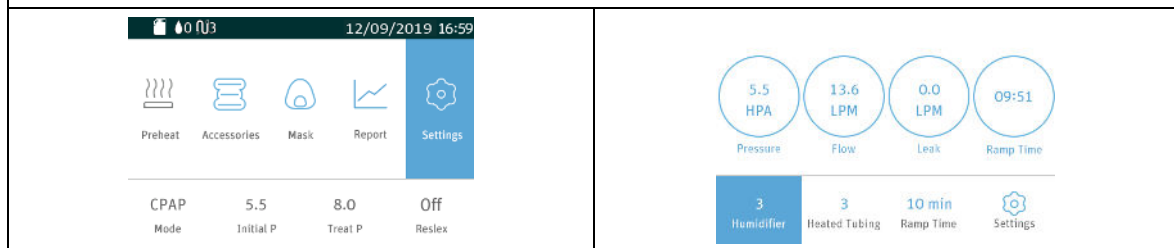
## 2. Setting up the Device





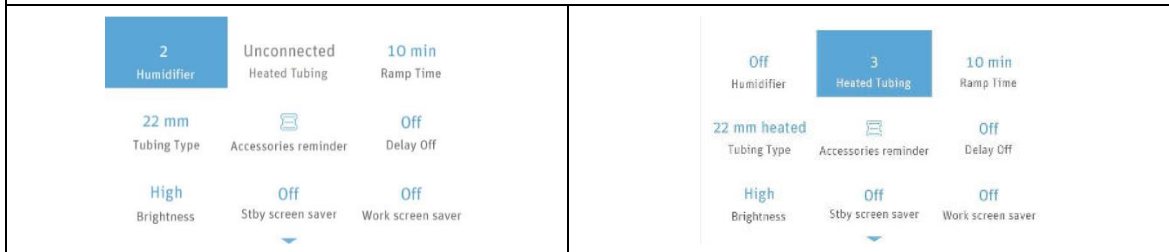
- A. Connect one end of the tube to the assembled mask.
- B. Connect one end of the tube to the air outlet of the device, as shown in the figure above. Or connect the heated tubing joint to the air outlet of the device, and then insert the power plug into the heated tubing port on the back of the device.
- C. Connect the power adapter to the DC inlet of the main device.
- D. Connect the power cord to the power adapter.
- E. Plug the power cord into the power outlet.


## 3. Adjusting Device Settings

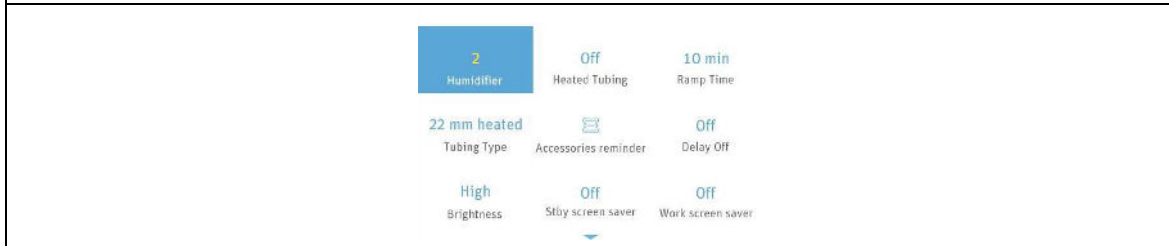
A. Connect the power cord and power adapter properly, The screen displays the Main Interface shown in the figure left below, press **the Start / Stop Button** , and the device will start delivering air, the screen displays the Main Interface shown the figure right below.



B. When the cursor is on the icon , press **the Knob**  to access the Initial Setup Interface of the Patient Menu, as shown in the figure left below. The **Heated Tubing** option can only be adjusted when the device is connected to the heated tube, as shown in the figure right below.




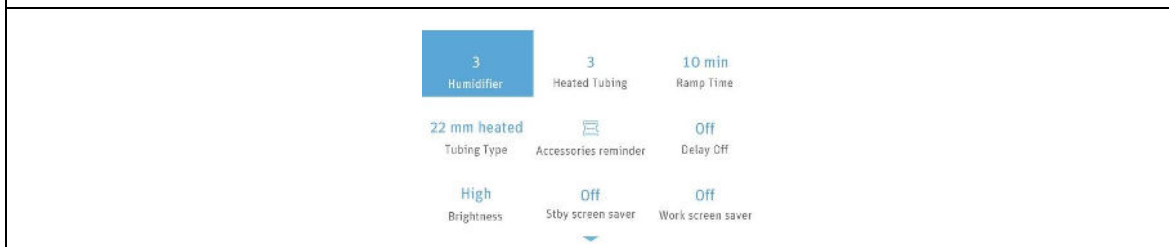
C. Press **the Knob**  to select option which is going to be changed, as shown in the figure below.





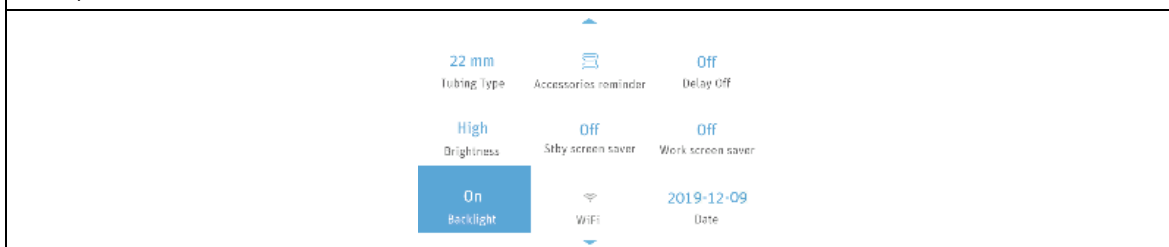
D. Turn **the Knob**  to adjust an option.




E. Press **the Knob**  to confirm your adjustment to the option.



F. The cursor switches among the options as you turn **the Knob** , as shown in the figure below. Refer to steps C-E to adjust other options. Press **the Home Button**  to exit the Setup Interface.



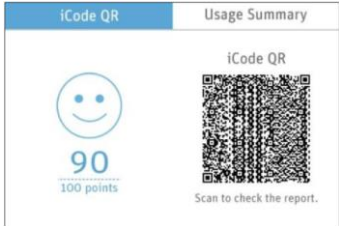
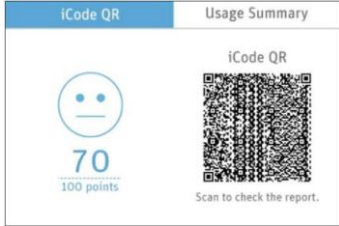
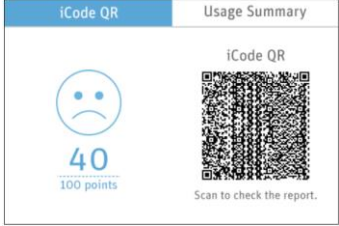
## 4. Starting Treatment

After you finish adjusting device settings and return to the Main Interface, press **the Start / Stop Button**  to start treatment.

## 5. Usage score of Luna® G3 BPAP 25A

After each treatment, Luna® G3 BPAP 25A will provide a Usage Score (for CPAP and AutoCPAP mode only) and Usage Summary on the screen for the user to get a quick understanding of the usage for the recent one day (from 12:00 noon of previous day to 12:00 noon of the day).

The total score is 100 points, with the below breakdown:

Total score	Status	Display
$\geq 75$	Good	
55 - 74	Normal	
$\leq 54$	Need to improve	

### Used time (60 points)

The more time of treatment the user received, the more points the user scores. If the device is used for 30 minutes, get 5 points, or for 1 hour, get 10 points. The user can score up to 60 points.

Used time (min)	$\geq 360$	330-359	300-329	...	60-89	30-59	0-29
Score	60	55	50	...	10	5	0

### Mask seal (20 points)

If the user is wearing the mask well without extra leakage during treatment, then the user can score up to 20 points. If the average leakage is less than 10 L/min, it scores 20 points, and 5 points will be deducted for every 5 L/min increase in leakage.

Leak (L/min)	<10	10-14	15-19	20-24	$\geq 25$
Score	20	15	10	5	0

### Respiratory events (10 points)

The fewer the respiratory events (apneas and hypopneas), the higher the score.

If AHI is less than 5, it scores 10 points, when  $AHI \geq 5$  and for every 5 increase in AHI, 2 points is deducted.

AHI	<5	5-9	10-14	15-19	20-24	$\geq 25$
Score	10	8	6	4	2	0